
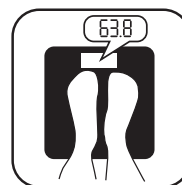
# OPERATING INSTRUCTION
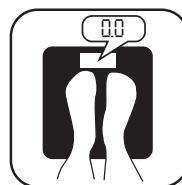
## GB Operating instructions for the scale

- Open battery cell, remove battery protection (plastic wrap). Select the kg/lb measurement unit using the switch located in the battery cell or in the middle of the scale's base.
- Set the scale on firm and level ground.
- Weight is measured after stepping on the scale.
- The measured weight will be displayed for 10 s. After that the scale turns off.
- Replace the battery if the indicator shows "Lo"
- If the indicator shows "Err", repeat the measurement

### COMMENTS:

1. The scale is designated **for household use only**.
2. Persons with a pacemaker should not use this scale.
3. To obtain correct readings weight should always be measured at the same time of day.
4. Do not step on the scale with wet feet to avoid slipping.
5. "Lo" on the display means low battery – in such a case replace the batteries.
6. "oL" on the display means that the allowed weight limit has been exceeded.
7. This product has been made using tempered glass. Store in a dry place not exposed to heat or humidity.
8. Shaking and throwing the scale may damage it. Glass breakage may result in cuts.
9. Use a moist cloth to clean the surface of the scale. Keep water from getting inside the scale. Do not use chemical agents.
10. **Do not place used batteries or electrical equipment in municipal waste containers.**

Take them to proper waste collection outlets or put them in specially marked containers.

### BASIC DATA:

- Load capacity: 5kg-150kg
- Gauge: 0.1 kg
- Working temperature: 10-35°C
- Batteries: 1 x 3V CR2032Li
- Display: LCD
- Scale units: kg/lb
- Low battery indicator: "Lo"
- Scale overload indicator: "oL"
- Error indicator: "Err"

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