Mom & Baby Mode

To weigh your baby, use Mom and Baby mode. First take your weight measurement. Then step off the scale, take the baby in your arms, step on the scale again. The app will show the baby's weight. If you want to add a child to the user list, follow the instructions.



Maintenance:

1. Clean with a soft, damp cloth.
2. Do not use chemical detergents for cleaning.

Important:

1. The scale is for home use only and measuring 14 parameters including BMI, fat, muscle, metabolism, hydration, bone mass.
2. It cannot be used for commercial services and other purposes.
3. Measurement of fat level is for reference only. If there are deviations from the norm, please contact your doctor or nutritionist.
4. Measurement range: 5-180 [kg].

Supported phone models:

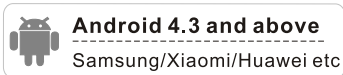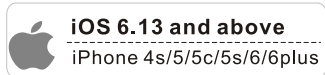
120UL/46.23

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Instrukcja

Electronic scale with analytical application

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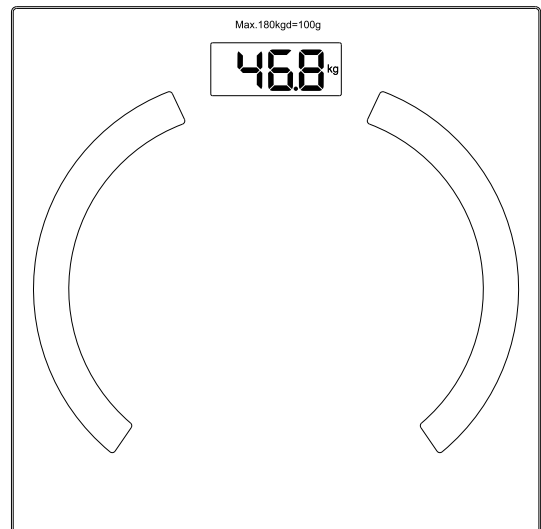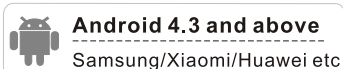
First launch:

Scan the QR code
(One code for iOS &
Android)

1. Starting up for the first time - after opening the package, install the batteries according to the instructions.
2. Download the AIFIT app (scan the QR code or download from Googleplay /Apple Store).
3. Install the app.
4. Open the app to register/create an account. Set up your data, turn on bluetooth to pair the scale with the app.
5. Get on the scale and take a measurement.

Measurements:

1. Get on the scale barefoot
2. Make sure your feet are touching the electrodes
3. The measurement takes 2-3 sec. (weight and BMI measurement). The result will show on the device paired with the scale.
4. After 5-8 sec. from the first measurement (point 3), the next results (fat, water, muscle, bone, etc.) will be displayed.

More features:

- You can add up to a maximum of 8 users (within one registered ID). The device recognizes a given user by parameters and historical measurements.
- When you open the application, the measurement data will synchronize with the device (scale).
- You can add your target (e.g. weight) to the device, follow the instructions.

ENGLISH

Package contents:

- Scale
- Manual
- 3 x AAA 1.5V batteries

CAUTION - before first use:

1. Do not use the scale if you have a pacemaker.
2. If you want to measure all parameters at the same time, you should step on the scale barefoot (without shoes or socks), otherwise you will only measure weight + BMI.
3. The scale may indicate incorrect measurement results when the measurement is taken for the first time or when the measurement location is changed. To get the correct result, turn off the device, turn it on again and take measurements again.
4. Measurements should be taken on a flat, hard surface. Measurements made on a soft surface such as carpet, mat, cardboard may not be correct.